

Original citation:

Dhaliwal, A., Zeino, Z., Tomkins, C., Cheung, M., Nwokolo, C., Smith, S., Harmston, C. and Arasaradnam, Ramesh P.. (2014) Utility of faecal calprotectin in inflammatory bowel disease (IBD) : what cut-offs should we apply? *Frontline Gastroenterology*, Volume 6 (Number 1). pp. 14-19.

Permanent WRAP url:

<http://wrap.warwick.ac.uk/65593>

Copyright and reuse:

The Warwick Research Archive Portal (WRAP) makes this work of researchers of the University of Warwick available open access under the following conditions.

This article is made available under the Creative Commons Attribution 4.0 International license (CC BY 4.0) and may be reused according to the conditions of the license. For more details see: <http://creativecommons.org/licenses/by/4.0/>

A note on versions:

The version presented in WRAP is the published version, or, version of record, and may be cited as it appears here.

For more information, please contact the WRAP Team at: publications@warwick.ac.uk

warwick**publications**wrap

highlight your research

<http://wrap.warwick.ac.uk>



RESEARCH

Utility of faecal calprotectin in inflammatory bowel disease (IBD): what cut-offs should we apply?

A Dhaliwal,¹ Z Zeino,¹ C Tomkins,³ M Cheung,¹ C Nwokolo,¹ S Smith,³ C Harmston,² R P Arasaradnam^{1,4}

¹Department of Gastroenterology, University Hospitals Coventry and Warwickshire, Coventry, UK

²Department of Colorectal Surgery, University Hospitals Coventry and Warwickshire, Coventry, UK

³Department of Biochemistry, University Hospitals Coventry and Warwickshire, Coventry, UK

⁴Clinical Sciences Research Institute, Coventry, UK

Correspondence to

Dr Ramesh Arasaradnam, Department of Gastroenterology, University Hospitals Coventry and Warwickshire, Clifford Bridge Road, Coventry, CV2 2DX, UK; R.Arasaradnam@warwick.ac.uk

Received 20 November 2013

Revised 8 February 2014

Accepted 26 February 2014

ABSTRACT

Background Faecal calprotectin (FC), a cytosolic protein released by neutrophils (S100 family) in response to inflammation, is a simple, non-invasive test that can be used to differentiate irritable bowel syndrome (IBS) with inflammatory bowel disease (IBD), where there can be considerable symptom overlap.

Aims and methods The aims of the study were (1) to be able to predict the ability of FC to exclude IBD and determine cut-offs when in remission, (2) to investigate the effects of time and temperature on stability of FC and (3) compare three ELISA kits to measure FC: Buhlmann, PhiCal v1 and PhiCal v2. A total of 311 patients with altered bowel habit were tested for FC; 144 with IBS, 148 with IBD and 19 with other organic causes.

Results Sensitivity and specificity of FC (with PhiCal v2 kit) to distinguish between functional disorder (IBS) and IBD using cut-off 50 µg/g were 88% and 78%, respectively, with a negative predictive value of 87%. Area under the receiver operating curve was 0.84 (CI 0.78 to 0.90). For those with IBD, FC values below 250 µg/g corresponded with remission of disease with a sensitivity and specificity of 90% and 76%, respectively. Area under the receiver operating curve was 0.93 (CI 0.89 to 0.97). FC was stable once extracted and frozen for up to 2.5 months. Pearson correlation was good between Buhlmann assay and PhiCal v2 ($r^2 = 0.95$).

Conclusions FC has up to 87% negative predictive value to exclude IBD, and cut-offs less than 250 µg/g had 90% sensitivity to determine remission in IBD. Once frozen, FC is stable and the ELISA monoclonal plates were broadly comparable.

INTRODUCTION

Faecal calprotectin (FC) is a calcium-binding heterodimer, which is abundant in the cytoplasm of neutrophils.¹ Inflammation causes neutrophil activation,

which results in calprotectin release proportionate to the degree of inflammation.² Gastrointestinal inflammation, as in active inflammatory bowel disease (IBD), releases calprotectin into body fluids including faeces.²⁻³ Manufacturers state that FC remains stable in faeces for up to 5 days³⁻⁴ and its use as a non-invasive biomarker in IBD management is gaining popularity. As calprotectin release is proportionate to the degree of inflammation, it may also have a role in monitoring disease activity.⁵ However, this requires local validation of the various available ELISA kits.

Different monoclonal ELISA kits are available for use in the UK, including the PhiCal (Immundiagnostik Calprotectin ELISA; R-Biopharm, Darmstadt, Germany) and Calprotectin ELISA (Buhlmann, Basel, Switzerland). The specification of the ELISA kits has changed over time as technology has evolved. A new version of the Immundiagnostik Calprotectin (PhiCal v2) ELISA was released in June 2011 replacing PhiCal v1. Previous studies, of each ELISA method, have shown correlation with IBD detection and monitoring of disease activity. Only one study to date has compared polyclonal and monoclonal ELISA for FC detection.⁶ However, no studies have compared the different sensitivities for disease detection or the reference ranges for disease activity of three monoclonal ELISAs. This is of importance when defining local reference ranges.

The aim of this study was to:

- ▶ Determine the ability of FC to exclude IBD and optimal cut-offs for IBD in remission
- ▶ To investigate the effects of time and temperature on stability of FC and
- ▶ To evaluate three different ELISA kits available commercially in the UK.

To cite: Dhaliwal A, Zeino Z, Tomkins C, et al. *Frontline Gastroenterology* Published Online First: [please include Day Month Year] doi:10.1136/fgastro-2013-100420

METHODS

This study was conducted at the University Hospitals Coventry and Warwickshire, which serves a population of approximately 500 000. Ethical approval was not required (as determined by the Local Research Ethics committee) as FC testing was already readily available in the UK.

A total of 311 FC samples were analysed prospectively between September 2010 and August 2012 in whom either a functional gastrointestinal (GI) disease (irritable bowel syndrome (IBS) predominant diarrhoea (IBS-D)) or IBD was suspected in a secondary care setting. As per protocol, all patients with IBS-D fulfilled the Rome 11 criteria (n=144) and details of patients with IBD were recorded (n=148); 57 with ulcerative colitis (UC) and 91 with Crohn's disease. There were 19 with other incidental conditions such as diverticular disease and microscopic colitis. All patients had screening for coeliac disease and anaemia.

Those with IBD had endoscopic/histological and/or radiological confirmation of disease and FC testing within 4 weeks. Demographic and clinical details are shown in table 1. Details of patients with IBD and FC levels are shown in table 2. Active disease was defined by presence of active inflammation on endoscopy, confirmed by histology; remission of disease was defined by quiescent disease on endoscopy and no active inflammation on histology.

To determine FC stability, stool samples from four random individuals were selected (one without symptoms and three with diarrhoea). These were received in the laboratory at room temperature within 6 h of voiding. Each was homogenised, then split into aliquots for immediate extraction to compare storage at room temperature and at 4°C for up to 3 days. Additionally, aliquots were also stored at 4°C, -20°C and -70°C and extracted at 1 week, 2 weeks and 6 weeks.

For comparison of ELISA kits, a random subset of nine individuals was analysed using three monoclonal ELISA plates; Buhlmann, PhiCal v1 and PhiCal v2. Three aliquots from each individual were analysed with the three ELISA plates in accordance with storage conditions described above. The updated PhiCal v2 (available only since June 2011) offers additional benefits with a shorter incubation period (30 min vs 1 h), room temperature incubation (versus 37°C shaking) and the use of a liquid stable conjugate.

Table 1 Demographics of patients with IBS and IBD

Groups (total; n=391)	IBS (n=144)	IBD active (n=88)	IBD remission (n=60)
M:F	1:1.9	1:1.8	1:2.9
CRP; mg/L (SD)	13 (4.6)	15 (17)	14 (43)
FC results; µg/g (SD)	34 (69)	674 (480)	92 (135)

CRP, c-reactive protein; FC, faecal calprotectin; IBD, inflammatory bowel disease; IBS, irritable bowel syndrome.

Table 2 Details of FC results in patients with IBD (UC and Crohn's) when active and in remission

Total IBD N=148	Crohn's (active) N=57	Crohn's (remission) N=34	UC (active) N=31	UC (remission) N=26
FC results; µg/g (SD)	646 (455)	96 (147)	649 (464)	108 (154)

FC, faecal calprotectin; IBD, inflammatory bowel disease; UC, ulcerative colitis.

This enabled a direct comparison between the old and the new kits by the same manufacturer.

FC sample analysis

For each sample, 100 mg of stool, using a 10 µL inoculation loop was weighed and dispensed into an analysis pot. The exact weight was recorded and 5 mL of extraction buffer added. Samples were vortexed for 30 min, in order to ensure dissolution, and then centrifuged using an Eppendorf centrifuge. The supernatant was removed for analysis by ELISA.

Data analysis and statistics

Data were compiled into Microsoft Excel for analysis and statistics was performed with SPSS V.17. Specificities, sensitivities, positive predictive values (PPVs) and negative predictive values (NPVs), for diagnosing IBD and when in remission, were calculated at different FC cut-off values. Non-parametric testing was applied to determine differences between groups.

Quality assessment

To ensure quality control across the ELISA plates, internal quality controls were run at the end of a plate and compared with the beginning, to minimise inter-observer variability. Additionally, samples were also sent for external quality assessment and local results were compared with published data.

RESULTS

To distinguish between IBS (non-organic disease) and IBD (organic disease), FC cut-off of 50 µg/g returned a sensitivity and specificity of 88% and 78%, respectively, with a PPV of 79% and a NPV of 87% (table 3). The area under the receiver operating curve (AUROC) was 0.84 (CI 0.78 to 0.90) (figure 1).

However, using cut-off 100 µg/g, sensitivity increased to 97% with a slight fall in specificity to 76% but with a PPV of 75% and a NPV of 97% (table 3). AUROC was 0.88 (CI 0.82 to 0.92) (figure 2).

For those with IBD, FC values below 250 µg/g correlated with remission of disease with sensitivity and specificity of 90% and 76%, respectively. AUROC was 0.93 (CI 0.89 to 0.97) (figure 3). Remission was

Table 3 Sensitivity and specificity at different cut-off ranges to distinguish IBS from IBD and IBD in remission compared with active disease

	Sensitivity (%)	Specificity (%)	PPV (%)	NPV (%)
IBS vs IBD at FC 50 µg/g	87.5	77.7	79.2	86.5
IBS vs IBD at FC 100 µg/g	97.2	76.3	75.3	97.4
IBD remission vs active at FC 50 µg/g	55.2	98.9	97.0	77.4
IBD remission vs active at FC 100 µg/g	72.4	95.6	91.3	84.3
IBD remission vs active at FC 250 µg/g	89.7	75.6	70.3	91.9

FC, faecal calprotectin; IBD, inflammatory bowel disease; IBS, irritable bowel syndrome; NPV, negative predictive value; PPV, positive predictive value.

confirmed by quiescent disease on endoscopy and no active inflammation on histology.

FC levels fell significantly when stored at room temperature compared with storage at 4°C, especially levels of FC <100 µg/g. If frozen at -20°C, then FC remains stable even prior to extraction into buffer solution for at least 6 weeks (figure 4). However once extracted into buffer solution, it remains stable for up to 2.5 months at -20°C.

Figure 5, shows FC results comparing the three ELISA monoclonal plates which were broadly comparable. There was good correlation between the Buhlmann and the PhiCal v2 results, with $r^2=0.95$ (for FC samples >250 µg/g) and $r^2=0.72$ (for FC samples <250 µg/g). The PhiCal v2 has a higher upper detection limit (840 µg/g) and greater linearity compared with the Buhlmann ELISA (600 µg/g); data not shown. External quality assurance for our results was comparable with the regional reference centre (figure not shown). The cost of ELISA kits are broadly comparable with PhiCal v2 costing ~£25/test.

DISCUSSION

A recent systematic review⁷ and meta-analysis⁴ have confirmed the utility of FC as a screening tool to identify patients with gut inflammation. Interestingly

in the meta-analysis only 30% of adults had confirmatory endoscopic or histological data to confirm inflammation. Our results demonstrate the utility of FC to distinguish non-organic (IBS) from organic (IBD) disease and provide insight into optimal cut-offs to determine when IBD is in remission.

Applying a cut-off level of 50 µg/g demonstrates acceptable sensitivity and specificity to distinguish between IBS and IBD. This does improve further when a cut-off of 100 µg/g is applied. The systemic review by Waugh *et al* which informed the National Institute of Health and Care Excellence diagnostics assessment group⁸ suggests most of the evidence is based on a 50 µg/g cut-off which enables reduction in the number of false negatives and is cost-effective.

Previous smaller studies have shown that FC concentrations correlate with endoscopic findings.² Thus, the ability to quantify FC with different severity levels of inflammation enables monitoring using FC to determine treatment response or failure in subjects

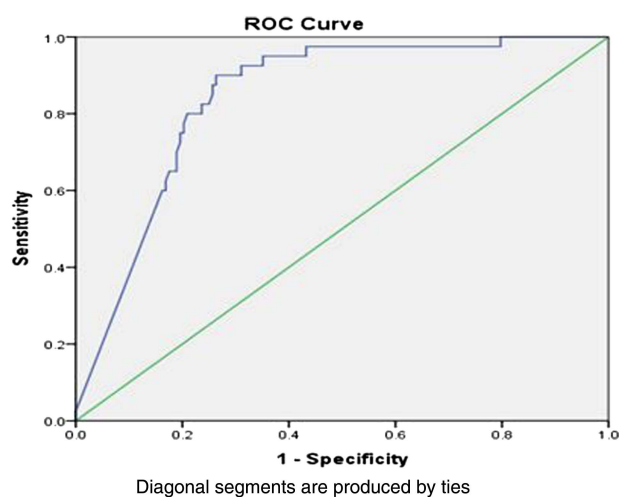


Figure 1 Area under the receiver operating curve showing 0.84 prediction for irritable bowel syndrome versus inflammatory bowel disease using cut-off at 50 µg/g.

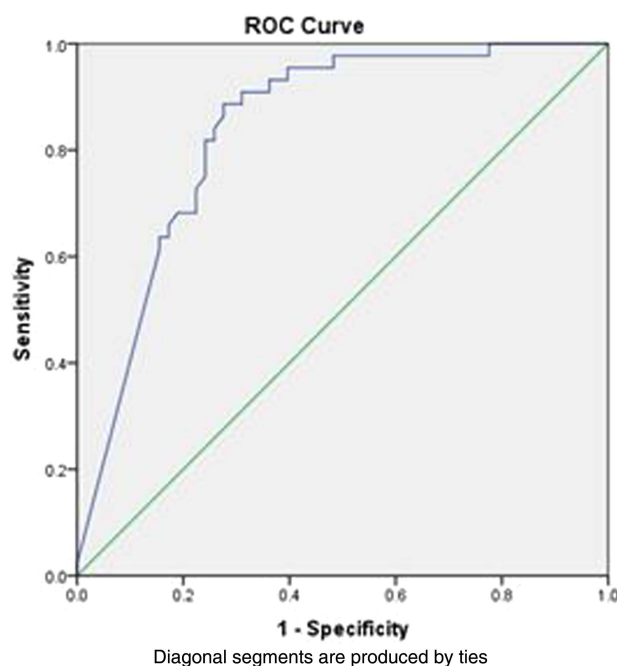


Figure 2 Area under the receiver operating curve showing 0.88 prediction for irritable bowel syndrome versus inflammatory bowel disease using cut-off at 100 µg/g.

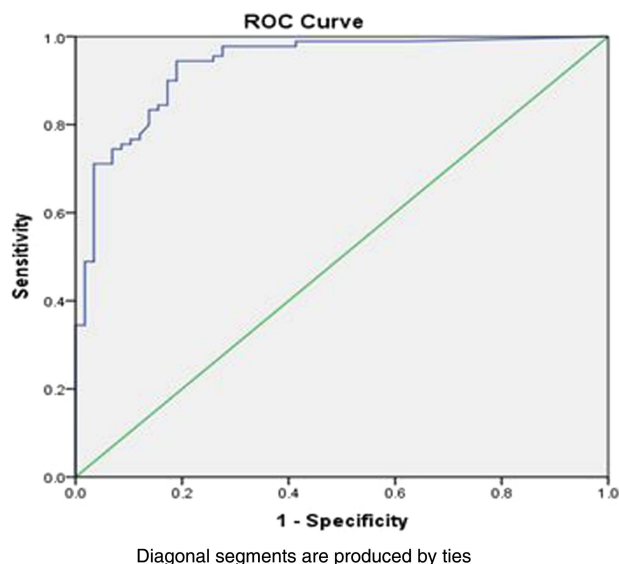


Figure 3 Area under the receiver operating curve (AUROC) for inflammatory bowel disease, active versus remission ROC was 0.93 (CI 0.89 to 0.97).

with active disease who are initiated on new therapy. This may ultimately reduce the need for repeated endoscopic assessments in active colitis. Damms and Bischoff,⁹ in 18 patients, have shown that the mean FC level in those with active IBD was 797. Mdel *et al*,¹⁰ in 25 patients with IBD, suggested a cut-off of 217 $\mu\text{g/g}$ to determine remission of IBD. In this study of 148 individuals, we have shown the optimal cut-off to be 250 $\mu\text{g/g}$ with 90% sensitivity and an AUROC of 0.93. FC levels for those with active Crohn's disease and ulcerative colitis were also broadly similar.

FC levels are influenced by storage conditions. Importantly however, once extracted and frozen, FC remains stable. Both manufacturers (Buhlmann and PhiCal v2) recommend a 'cut-off' FC level of 50 $\mu\text{g/g}$, above which all results are deemed positive. Our data also confirms the presence of improved linearity at higher values. This has specific relevance

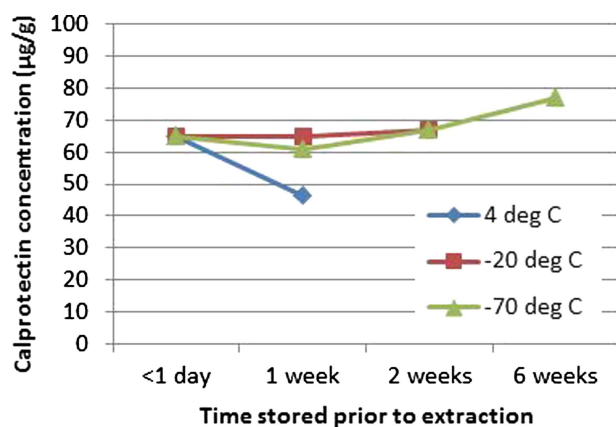


Figure 4 Faecal calprotectin levels based on temperature against time prior to extraction into buffer solution.

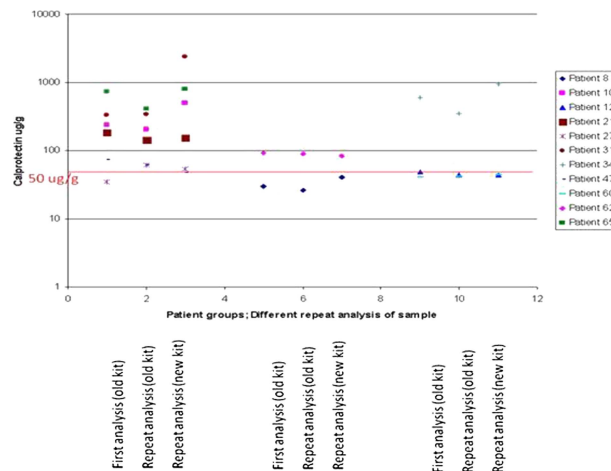


Figure 5 Comparison of faecal calprotectin values with three ELISA kits (Column 1—Buhlmann, Column 2—PhiCal v1 and Column 3—PhiCal v2).

when monitoring disease activity and treatment response in subjects known to have IBD. On a technical note, the PhiCal v2 ELISA kit has a higher upper limit of detection thus reducing the number of dilution steps for high FC values. This would be an important consideration in monitoring IBD disease activity but less so for distinguishing between IBS and organic disease/IBD. Compared with the older PhiCal v1 (no longer available, superseded by PhiCal v2), it has a shorter incubation period and easier application.

It is worth noting that for the same patient groups the Buhlmann ELISA kit consistently produced higher numerical (~20%) FC results (figure not shown). This has led some departments to consider revising the cut-off value to, for example, 100 $\mu\text{g/g}$. Intriguingly, both ELISA kits use the same monoclonal antibody, thus the differences must lie in the buffer solution at the time of extraction, reinforcing our views that this appears to be the critical step when analysing FC.

Although all FC testing (even with different ELISA kits) was performed manually, our findings were comparable with other reported series^{1–6} supporting the validity of our results. Moreover, internal and external quality analyses of our technique were acceptable given the low sample numbers and manual pipetting techniques—within the initial validation cohort of 62 individuals. Subsequent automation for improved precision enabled throughput of larger sample numbers with improved coefficient of variance within our laboratory.

C-reactive protein (CRP) and erythrocyte sedimentation rate (ESR) are laboratory investigations which are often employed by clinicians in their determination of probability of patients having an underlying organic bowel disease as opposed to a functional disorder. Using a low CRP cut-off is reported to have improved sensitivity, but poorer specificity, which would lead clinicians to unnecessarily investigate, if

used in isolation.^{11 12} ESR has been shown, in other studies, to offer little benefit to the evaluation process for those subjects fulfilling clinical criteria suggestive of IBS.¹² Tibble *et al* reported that abnormal FC had a significantly greater OR using FC to identify IBD ($p < 0.0001$) compared with an elevation in CRP and ESR.¹³ Our findings support this earlier result by showing greater sensitivity and predictive ability with FC than with ESR and/or CRP to distinguish between functional and organic disorders. This in turn, aids in monitoring of active exacerbations and remission of IBD.

FC is more cost-effective than non-FC strategies in a diagnostic pathway based on a recent systematic evaluation by the Centre for Economic Based Practice 2010 data.¹⁴ Although any quality-adjusted life-year gains were likely to be small due to of the low prevalence of IBD and the high sensitivities of all of the tests, considerable savings could still accrue.⁷ Additionally, a recent study of 3639 patients has suggested that the estimated demand for colonoscopies is reduced by 50% with the 50 µg/g cut-off and by 67% with the 100 µg/g cut-off, equating to €1.57 million and €2.13 million, respectively.¹⁵

Overall our data affirms that the use of FC is beneficial in distinguishing between functional GI conditions (IBS) and organic disease (IBD). In those with IBD, a 250 µg/g cut-off aids in determining clinical disease activity. FC testing like β natriuretic polypeptide used in the diagnosis of cardiac failure, is a useful, relatively economical test benefiting primary and secondary care. It helps to sift out those with organic pathology requiring further investigations, and in those with IBD, aids in disease monitoring, avoidance of invasive procedures and reduction of time off work.⁷ The study population was carefully chosen to allow us to validate symptoms and assess all subjects with endoscopy. Further studies are now needed to assess the use in a primary care population.

What is already known on this topic

- ▶ Faecal calprotectin (FC) has a high negative predictive value to distinguish between organic and non-organic disease (functional).

What this study adds

- ▶ FC levels fall at room temperature if stored inappropriately ELISA kits to - measure FC are broadly similar; the extraction process is critical.
- ▶ For those with inflammatory bowel disease (IBD), values of >250 µg/g indicate disease activity.
- ▶ FC values for active IBD are not lower in Crohn's disease compared with ulcerative colitis.

How might it impact on clinical practice in the foreseeable future

- ▶ FC testing which is non-invasive, can guide clinicians in managing patients with inflammatory bowel disease.

Contributors All authors contributed to this work. RPA conceived and designed the study. AD, ZZ, MC, CH contributed to writing of the manuscript and data analysis. CT performed the data analysis of faecal calprotectin. RPA, CH, SS and CN revised the paper for intellectual content. All authors discussed the results and implications and commented on the manuscript at all stages.

Competing interests None.

Patient consent Obtained.

Provenance and peer review Not commissioned; externally peer reviewed.

Data sharing statement Further data is available on request from the corresponding author at R.Arasaradnam@warwick.ac.uk.

Open Access This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 3.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/3.0/>

REFERENCES

- Roseth AG, Fagerhol MK, Aadland E, *et al*. Assessment of the neutrophil dominating protein calprotectin in feces. A methodologic study. *Scand J Gastroenterol* 1992;27:793–8.
- Aomatsu T, Yoden A, Matsumoto K, *et al*. Faecal calprotectin is a useful marker for disease activity in pediatric patients with inflammatory bowel disease. *Dig Dis Sci* 2011;56:2372–7.
- Striz I, Trebichavský I. Calprotectin—a pleiotropic molecule in acute and chronic inflammation. *Physiol Res* 2004;53:245–53.
- van Rheenen PF, Van de Vijver E, Fidler V. Faecal calprotectin for screening of patients with suspected inflammatory bowel disease: diagnostic meta-analysis. *BMJ* 2010;341:c3369.
- Sipponen T, Kärkkäinen P, Savilahti E, *et al*. Correlation of faecal calprotectin and lactoferrin with an endoscopic score for Crohn's disease and histological findings. *Aliment Pharmacol Ther* 2008;28:1221–9.
- Malickova K, Janatkova I, Bortlik M, *et al*. Calprotectin levels in patients with idiopathic inflammatory bowel disease comparison of 2 commercial tests. *Epidemiol Mikrobiol Imunol* 2008;57:147–53.
- Waugh N, Royle P, Cummins E, *et al*. Faecal calprotectin testing for differentiating amongst inflammatory and non-inflammatory bowel diseases: systematic review and economic evaluation. *Health Tech Assessment* 2013;17:xv–xix, 1–211.
- NICE guidelines; Faecal calprotectin diagnostic tests for inflammatory diseases of the bowel; DAP 12; June 2013.
- Damms A, Bischoff SC. Validation and clinical significance of a new calprotectin rapid test for the diagnosis of gastrointestinal diseases. *Int J Colorectal Dis* 2008;23:985–92.
- Mdel GS, González R, Iglesias FE, *et al*. Diagnostic value of faecal calprotectin in predicting an abnormal colonoscopy. *Med Clin (Barc)* 2006;127:41–6. (abstract in English).
- Jellema P, van Tulder MW, van der Horst HE, *et al*. Inflammatory bowel disease: a systematic review on the value

- of diagnostic testing in primary care. *Colorectal Dis* 2011;13:239–54.
- 12 Langhorst J, Elsenbruch S, Koelzer J, *et al.* Noninvasive markers in the assessment of intestinal inflammation in inflammatory bowel diseases: performance of fecal lactoferrin, calprotectin, and PMN-elasticase, CRP, and clinical indices. *Am J Gastroenterol* 2008;103:162–9.
 - 13 Tibble JA, Sigthorsson G, Foster R, *et al.* Use of surrogate markers of inflammation and Rome criteria to distinguish organic from nonorganic intestinal disease. *Gastroenterology* 2002;123:450–60.
 - 14 CEBP economic report: Value of Calprotectin in screening our IBS 2010. http://www.alphalabs.co.uk/asp/db/documents/cep09026_Calprotectin%20Economic_Jan10.pdf
 - 15 Mindemark M, Larsson A. Ruling out IBD: estimation of the possible economic effects of pre-endoscopic screening with F-calprotectin. *Clin Biochem* 2012;45: 552–5.



Utility of faecal calprotectin in inflammatory bowel disease (IBD): what cut-offs should we apply?

A Dhaliwal, Z Zeino, C Tomkins, M Cheung, C Nwokolo, S Smith, C Harmston and R P Arasaradnam

Frontline Gastroenterol published online March 17, 2014

Updated information and services can be found at:

<http://fg.bmj.com/content/early/2014/03/17/flgastro-2013-100420>

These include:

References

This article cites 12 articles, 1 of which you can access for free at: <http://fg.bmj.com/content/early/2014/03/17/flgastro-2013-100420#BIBL>

Open Access

This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 3.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/3.0/>

Email alerting service

Receive free email alerts when new articles cite this article. Sign up in the box at the top right corner of the online article.

Topic Collections

Articles on similar topics can be found in the following collections

[Open access](#) (28)

Notes

To request permissions go to:

<http://group.bmj.com/group/rights-licensing/permissions>

To order reprints go to:

<http://journals.bmj.com/cgi/reprintform>

To subscribe to BMJ go to:

<http://group.bmj.com/subscribe/>